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DE

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## Over-bed Table RFM

UK




 Carefully read the complete user instructions before assembling, setting or using the Over-bed Table RFM for the first time.

### Intended Purpose

For people who are bed-bound due to illness or who, following an operation, cannot get out of bed too often to begin with, the over-bed table is ideally suited as an aid for eating and drinking and as a storage shelf. The castors on the bottom of the frame enable the individual positioning of the table next to the bed so that strenuous movements can be avoided. The setting of the height of the table serves for the reaching of objects or tableware, regardless of the size / height of the bed.

### Mounting / Operation

- The bed table RFM is delivered pre-assembled. Check the contents for completeness.
- Remove the lower frame and attach the wheels to it. Make sure that the wheels with brakes are located at the outer/rear ends of the frame.
- Put the upper stand upright with the table top and fix it to the lower stand with the black locking screw.
- There is a red clamping screw under the table top. To adjust the inclination of the table top, loosen it. This loosens the tothing. Then adjust the table top to the desired inclination and tighten the clamping screw again. Make sure that the tothing closes again.
- You can adjust the height of the bed table using the knurled screw on the frame.
- For a safe stand, lock the two brakes on the rear wheels.
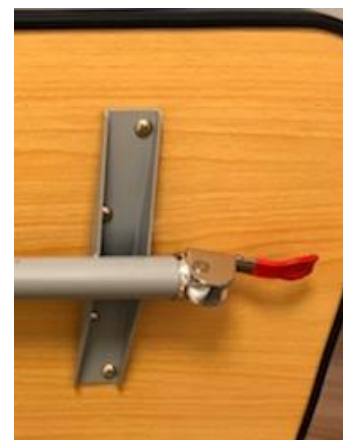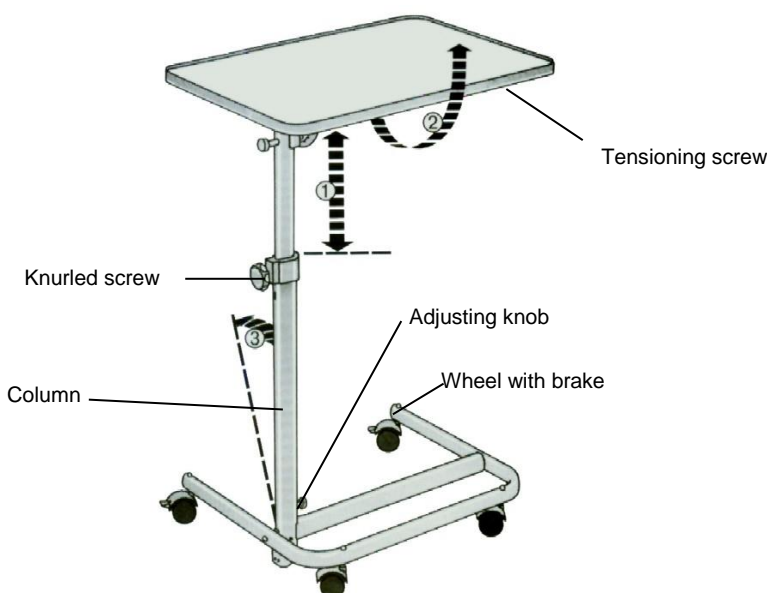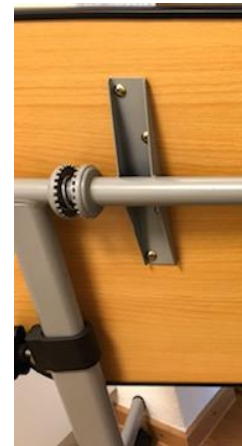
## Safety Instructions

- The Over-bed Table RFM is intended for domestic use and before every use it must be checked that the knurled screw and the locking screws are securely tightened!
- The correct operation of the brakes must be checked regularly.
- **Never use a defective over-bed table!** Please contact your specialist dealer or supplier if the material is missing or has been altered.
- Do not support your body weight on the table top and don't pull yourself up or into a sitting position using the Over-bed Table RFM. The bed table can tip over and is only intended as a shelf for objects! The maximum load-bearing capacity is 18 kg!

## Cleaning / Multiple Use

Regularly clean the Over-bed Table RFM with a mild soap solution and clean water. Use a soft cloth for drying. Use a commercially available disinfectant to disinfect the bed table.

The RFM bed table is intended for daily use. It can be used by several persons. Before the table is used by another person or becomes the property of another person, it must be thoroughly cleaned and disinfected. Carry this out as described above.

## Maintenance

The Over-bed Table RFM is maintenance-free. If you notice that brakes on the castors weaken, push buttons do not fully lock into place or the screws cannot be fully tightened, contact your specialist dealer!

Do not use a defective bed table!

## Storage

Protect the Over-bed Table RFM against damp and store it in a safe place to protect the frame, the table and the castors against damage!

## Scope of Supply

Over-bed Table RFM

User Instructions

Outer packaging / Polybag

## Technical Data / Dimensions and Weight

Width:	61 cm
Depth:	45 cm
Table top:	62 x 41 cm
Weight:	7.5 kg
Size folded:	57 x 85 cm
Height adjustment:	72 – 94 cm
Maximum load:	18 kg
Material:	Painted steel (frame), plastic / PA (castors), plastic / coating (table top / safety edge)

## Guarantee

The guarantee is for 24 months and relates only to product defects which can be verifiably traced back to manufacturing faults or material defects.

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