Pedaltrainer RFM[®] Deluxe

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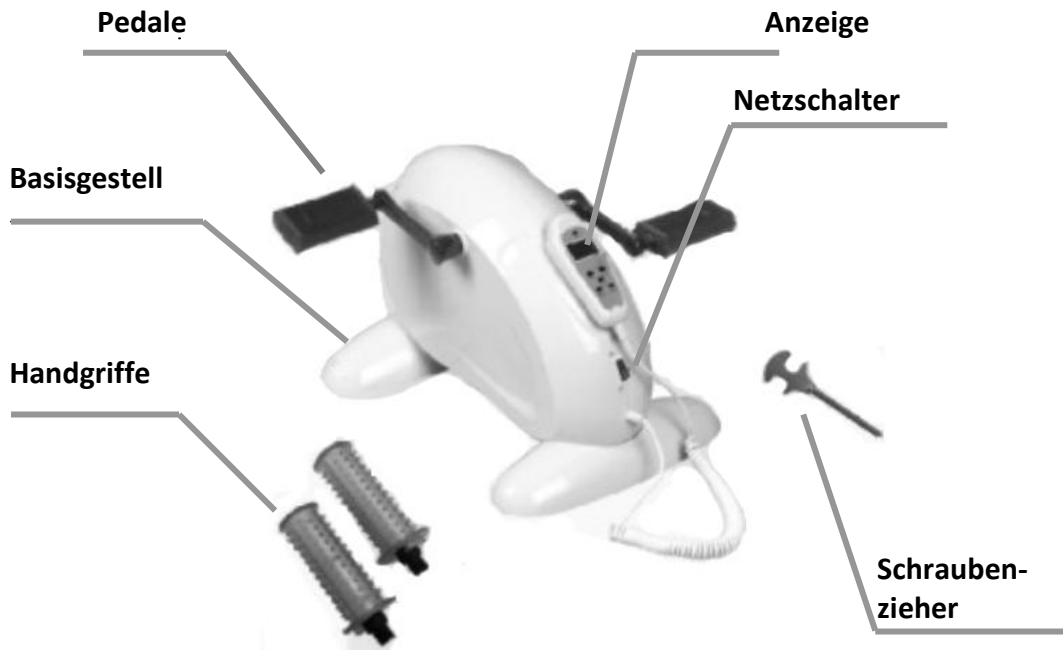
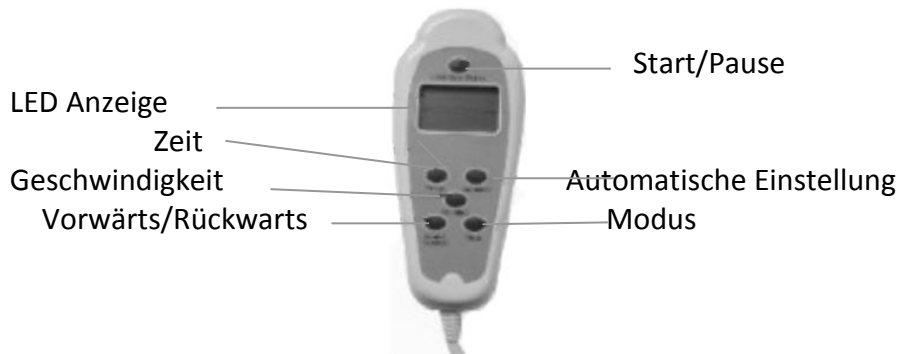
Pedaltrainer RFM[®] Deluxe

Important safety instructions

If you use an electrical device, you should always heed the basic precautionary measures including the following:

Read all instructions before use.

Danger – to reduce the risk of electrical shock:

Unplug the device immediately after use and before you clean it.

Warning– to reduce the risk of burns, fire, electrical shock or personal injury:

1. Never leave the device unsupervised when it is plugged in. Unplug it when not in use, and when you attach or remove parts.
2. Do not use under a blanket or pillow. Excessive heat can cause fire, electrical shock or personal injury.
3. Watch the device carefully when used in the proximity of children, sick or handicapped persons.
4. Only use the device for its intended purpose as described in these operating instructions. Do not use any accessories that are not recommended by the manufacturer.
5. Do not use the device if the power cord or plug is damaged. If it does not function properly, was dropped or damaged or fell in water, return it to customer service for inspection and repair.
6. Do not carry this device by the cord and to not use the cord as a handle.
7. Keep the cord away from hot surfaces.

8. Do not use the device if the vents are blocked. Keep the vents free of lint, hair, etc.
9. Do not let objects fall into the openings and do not stick anything in them.
10. Do not use the device outdoors.
11. Do not use the device near aerosol spray cans or areas where oxygen is in use.
12. To unplug the device, first switch all switches to *off* and then pull the plug from the outlet.
13. Keep children away from the extended pedals.

Features

This product is only intended for household use.

Normal operation of this device:

- a. Speed control: 4 different speeds 30/40/50/60 R/min
- b. Time settings: 5 time settings 3/6/9/12/15 min
- c. Forward and backward function: wait 2 seconds before changing direction
- d. Display function for overall speed reached
- e. Calorie counter
- f. Automatic run function: when the button is pressed in normal position, the fun function begins automatically
- g. Signal emitter: if this device is incorrectly operated, the motor will shut down automatically and the buzzer will sound an alarm; press once on the on/off switch and you can re-start the device – if you do not release the button after one minute the device will automatically shut down
- h. No stress is placed on the hips and knees

Warning

1. When not in use turn the switch to *off*.
2. Never submerge the device in water and never wipe with a wet cloth.
3. Do not clean it with a chemical cleaning agent.
4. Do not use it if it is not operating properly.

Operator's instructions

1. Make sure that the voltage is set to 220V
2. Attach pedal the pedals <Fig. 1>
3. Attach the handles <Fig. 2>
4. Plug the device in and switch it on; the contrast signal consisting of 4 different 8 digit numbers should appear on the display; – then set the functions, e.g. speed 30 R/min; time: 9 min; direction; display RPM – then press the start/pause button and you can begin

Fig. 1

Fig.2

Button function

1. Start/pause button: press once to start –press again and hold for 2 seconds to stop
2. Speed control button: press the speed button – 30/40/50/60 R/min can be selected
3. Forward and backward function button: Change direction by pressing once
4. Time setting button: 3/6/9/12/15 minutes – 5 different time units can be selected
5. Automatic button: The motor will run according to the selected program
6. Mode button: Display RPM, calories, meters
7. Signal emitter: If the device is operated improperly, the motor will stop automatically and the buzzer will sound an alarm; press once on the on/off switch, then you can re-start the device- if you do not release the button after one minute's time, the device will shut off automatically.

Instructions for Accessories

1. Shut off the device first, if you are attaching or replacing accessories
2. Assemble the pedals according to figure 3
3. Assemble the handgrips according to figure 4

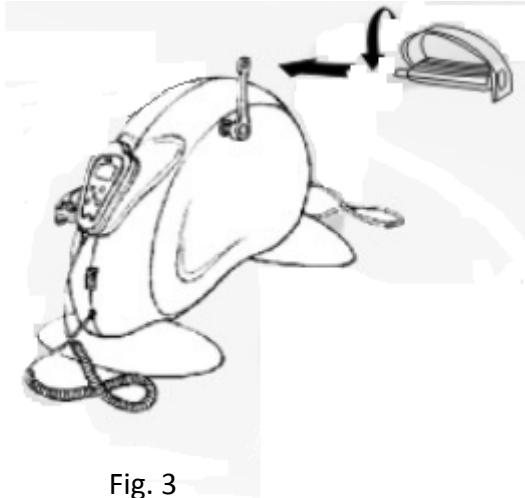
Fig. 3

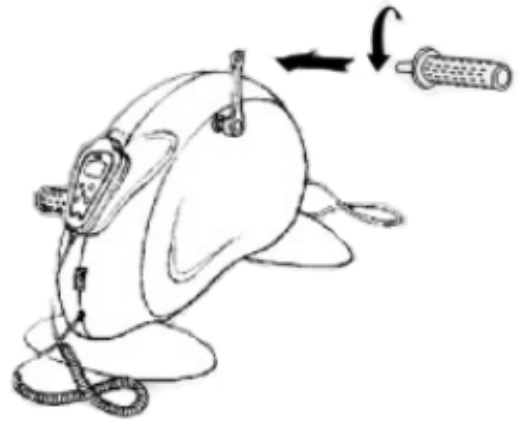
Fig. 4

Maintenance

This device is for household use only. There are no parts requiring maintenance inside. All maintenance should be performed by a recognised service provider. If cleaning is necessary, unplug the device and clean only the external surface with a cloth. Never place the device in water.


Make sure that the cord is not strained or twisted. Do not wrap the cord around the device.

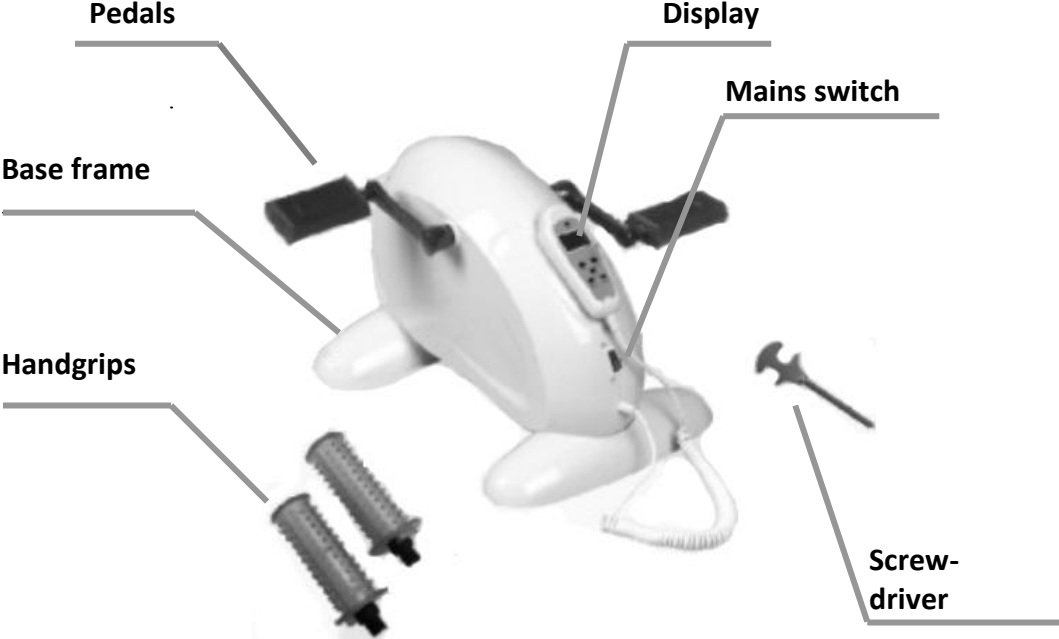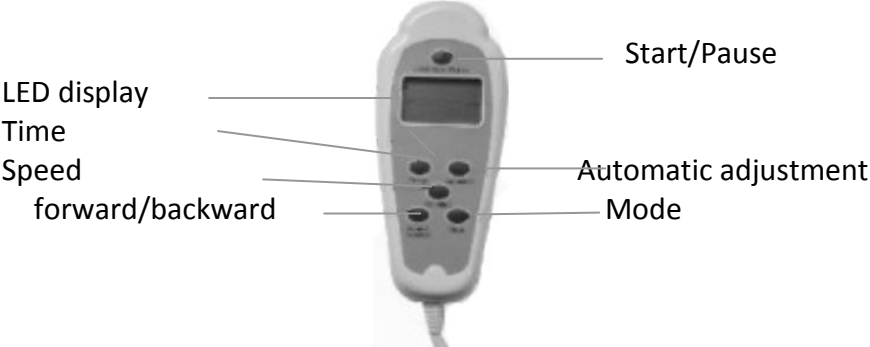
When not in use, the device should remain unplugged and stored in a safe, dry place— out of the reach of children. Do not wrap the power cord around the device.

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