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### **SAFETY PRECAUTIONS**

1. People with the following conditions should consult their physician before using the product:
  - Who uses a medical electronic device, which is embedded in the body, such as a pacemaker.
  - Who suffers from heart disease.
  - Who suffers from backbone injury or disease.
  - Who is suffering from bone fragility, especially on the backbone.
  - After surgery.
  - Who is pregnant or in a period of monthly illness.
  - Who suffers from malignant tumors.
  - Who is suffering from osteoporosis.
  - Who is suffering from blood circulation problems or thrombosis.
  - Who is suffering from skin disease or infection.
  - Who is undergoing medical treatment or feels some physical abnormality previously or is being treated now.
  - Who has some abnormality or deformity of limbs or backbone.
  - Who has injuries, diseases due to illness or accident.
2. Do not let water enter the product or be splashed near the surface of the product.
3. Water contact with electricity is dangerous. Do not leave the product in a humid environment (such as in the bathroom or shower room or swimming pool).
4. Before connecting the product to the power supply, check that the product voltage and supply voltage are the same.
5. If the plug, power cord or product itself is damaged, do not use the product.
6. In order to prevent the risk of electric shock, do not open the product outer cover.
7. Product is not intended for use by the physically or mentally impaired (including children), unless there is supervision or guidance by the people responsible for their safety when they use the product.

8. If the product or the power cord is damaged, do not try to repair the product. This must be handled by OTO Bodycare full-time specialist / technician otherwise the product guarantee will become void.

9. Children should not use this product. They should be under supervision to ensure that they will not treat this product as a toy.

10. If an unusual sound or noise is being heard from the product, turn off the product and immediately disconnect the power connection. Please contact an authorized OTO Service center immediately.

11. Do not use the product if cover is damaged or if there is a mechanical or electronic failure. Contact the authorized OTO Service center.

12. Only use the product in accordance with the instructions in this user manual.

13. Do not use this product in high temperature conditions or beside hot objects such as heaters, stoves.


14. Do not apply medical or cosmetic creams on the product or on your foot before using this product.

15. Do not destroy, bend stretch, twist or cut the power cord to avoid electric shock.

16. People who are sensitive to heat should not use the heat function of product.

17. Recommended usage time for an individual is 15 minutes at a stretch.

### **PRODUCT DESCRIPTION**

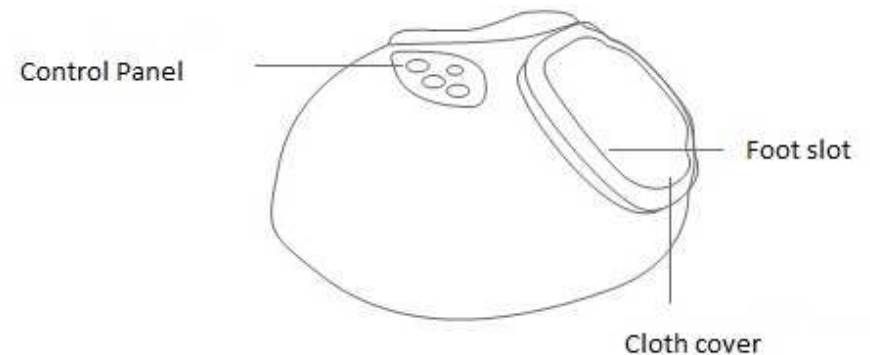**CONTROL PANEL**

**A. MODE**

“AUTO” Mode : A Combination of Kneading and Air Massage for Feet


“MASSAGE” Mode : Kneading Massage for Soles

“AIR” Mode : Full Air Massage for Feet

**B. INTENSITY levels**

**C. HEAT Therapy**

**D. STOP**

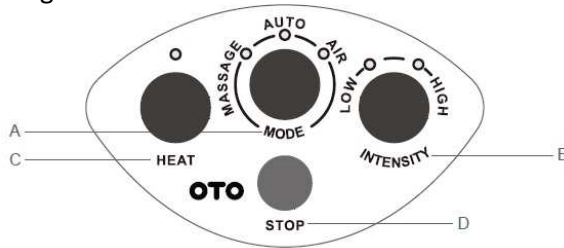**Key functions:**

Button	Function
MODE	Press to turn ON with default “AUTO” mode of low intensity Press again to turn on “MASSAGE” mode Press one more to turn on “AIR” mode Press again to stop all massages
STOP	Press to stop the massage any time
INTENSITY	2 Levels of intensity: HIGH / LOW (Applicable to “AUTO” and “AIR” modes only)
HEAT	Heat therapy of the soles

Auto Timer: 15 minutes. Power will turn Off automatically after 15 minutes to prevent over heating.

Always unplug the AC power cord after use.

**TROUBLE-SHOOTING**

This page summarizes the most common problems you could encounter with the product :

Problem	Possible Cause	Suggestion
The product does not work	The product is not plugged in properly.	Ensure the voltage on the product corresponds to the local mains voltage and then plug in again.
	The product has just session and has automatically switched to standby mode.	If you want another session, press the MODE button again.
	You may have used the product continuously for 30 minutes or more where the over heat protection has been activated.	Unplug the product and let it cool down for 60 minutes before you use it again.
	An excessive load has been applied on the product.	Reduce the excessive load and start again.
The product is abnormally warm	The product is likely to be overused.	Unplug the product and let it cool down before using it again.
Noise can be heard from the product	Noise is made by the mechanical parts of the product.	This is normal and does not mean that the product is faulty.

M1501 / KD1501

## **CARE AND MAINTENANCE**

### **Cleaning**

Do not use steel wool, abrasive cleaning agents or corrosive liquids (such as gasoline or acetone) to clean the product.

1. Unplug the plug from the power outlet.
2. Wipe with a dry cloth to clean the product.
3. Do not use a damp cloth or detergent

### **Storage**

1. Unplug the plug from the power outlet.
2. Store the product in a cool dry place away from sunlight.
3. Do not put a heavy weight load on the product.

If the cover is damaged or there is any mechanical failure, please stop using the product.

The European directive 2012/19/EC about electronic and electric waste, requires that you can't throw away defective domestic appliance with common waste. Used device has to be collect separately to optimize the recuperation rate and the recycling of the materials to reduce the impact on the health and the environment.

*This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capacities, or by persons lacking in experience or knowledge, unless they are supervised or instructed on the use of this apparatus by a person responsible for their safety. Children should be supervised to ensure they do not play with this appliance.*

### **Precision and precautions about the battery use**

- Throwing away a used battery in the environment or with household waste pollutes and avoids the recuperation of reusable materials.
- It is important to limit the battery consumption and to respect the following rules:
  - Favor alkaline battery (last longer than saline battery) and when it's possible use rechargeable battery.
  - Leave used batteries and accumulators in adapted containers that you can find in some retailers. In that way, metals don't pollute the environment and the health.
- Battery has to be inserted by respecting the polarity show on the apparel and the battery. An incorrect positioning can damage the apparel, cause leaking and even provoke fire and the explosion of the battery.
- To assure a good functioning, battery has to be in a good condition. If there is a problem, put new battery.

- Never try to recharge battery not made for this. It could leak, warm up, cause fire or explode.
- Replace all the battery at the same time. Never mix up alkaline battery with saline battery.
- Used battery has to be removed.
- Also, remove the battery if you're not using the apparel for a long time.
- Never try to short-circuit the terminal of a battery.
- Never throw away the battery in fire because they can explode.
- Take out the accumulator from the apparel before reload it.
- We recommend that an adult watch if a child change the battery to respect those rules.
- If you swallow a battery, go immediately to your doctor or the hospital.
- For the best ongoing performance, we recommend using only fresh, high-quality, alkaline batteries.

**Importé par / Importiert von / Geïmporteerd door / Imported by :  
CAPTELEC, 59170 Croix - FRANCE**