

MM384/ V7200 MASSAGING AND HEATING BED CUSHION WITH LED LAMP

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.



This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.

Children must be supervised to ensure that they do not play with the device.

SAFETY INSTRUCTIONS

1. This appliance is intended for private and non-commercial use.
2. Use the appliance only for its intended use described in this manual.
3. Before connecting the appliance to the mains, make sure that the power cord and plug are in good condition and that the voltage indicated on the appliance matches that of your electrical installation.
4. Only use with the supplied adapter.
5. Never leave the appliance unattended when plugged in.
6. **DO NOT USE** the bed cushion if it is folded.
7. Check that the cable is not compressed by the device or by any heavy objects.
8. Do not wrap the cable around the bed cushion.
9. After use, unplug the adapter from the power outlet. Do not pull on the cord, preferably grasp the plug.
10. Before each use, carefully check the condition of the adapter, cables, function remote control and device.
11. Do not use the device when the room temperature exceeds 40° C.
12. Keep it away from heat sources, direct sunlight, humidity, water or other liquids.
13. Do not allow the device to come into contact with pointed or sharp objects.
14. Never immerse the appliance in water or any other liquid.
15. **DO NOT USE** the device if your body is wet or if you have transpired.
16. **DO NOT USE** in damp rooms or near liquids, eg in a tub or in a sink filled with water, under a shower or in a pool.
17. **DO NOT USE** the device if one of its elements are damaged, if they are wet, if water or liquids penetrate inside the device or if it does not function properly.
18. If the device is very damaged or some of its parts are made bare or torn, do not use it. Turn it off immediately and disconnect the plug. Using the device in such cases could cause injury resulting from electric shock.
19. Do not repair, disassemble or modify the device and its elements.
20. If the product is damaged, it must be repaired by a qualified technician, in order to avoid all risks.
21. Do not subject it to strong impacts and do not throw it.

HEALTH PRECAUTION




Please consult your doctor prior to use if any of the following describes your physical condition

- Pregnant
- Suffer from heart problems, especially if you wear a pacemaker
- Benign or malignant tumors
- If you suffer from the following diseases or illnesses : epilepsy, muscle or joint problems, diabetes (especially if you take insulin), vascular problems, circulatory problems, varicose veins, open wounds, bruises, crack skin and inflammation of veins.
- If you have unexplained pains, if you follow a medical treatment, if you use medical devices , if you have or had serious illness not specified in this instruction manual and in case of ailments requiring special precautions.

- This device is not designed for medical use. sThanks to its two functions: massage and heating effect, it allows wellness (massage) and sweating action (heating effect).
- Do not use while sleeping. If you feel signs of drowsiness, turn it off.
- Do not use 30 minutes before or after having a meal.
- Persons who are insensitive to heat should only use massage function.

PRODUCT DESCRIPTION

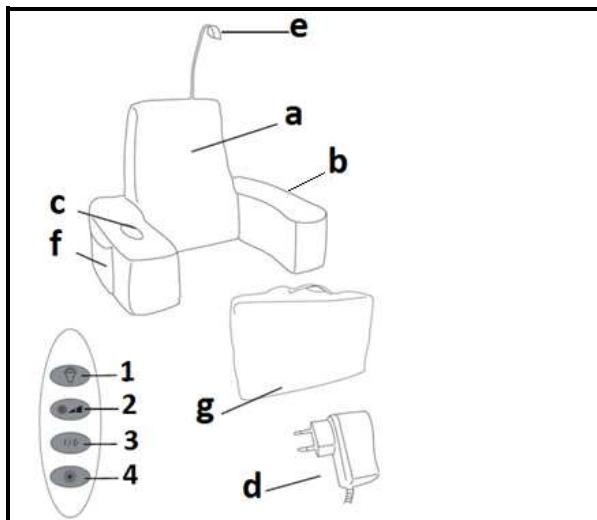

- a. Backrest
- b. Folding armrests
- c. Fuction remote control
 1. Light function
 2. Massage intensity
 3. I/ O ON/ OFF massage
 4. Heat function
- d. Adapter
- e. LED lamp with swivel arm
- f. Side pockets
- g. Slipcover

TECHNICAL SPECIFICATIONS

Adapter :

- Input: 100-240V ~ 50/ 60Hz 0.5A
- Output: 12V --- 0.8A

Massaging and heating bed cushion:

- 12V --- 0.8A 9.6W

OPERATION

1. Correctly mount the unit.
2. Place it where you want to use it : bed or sofa.
3. Sit down.
4. Connect the adapter cable to the unit.
5. Familiarize yourself with the different functions of the function remote control: Light, massage or heating function.

OPERATING THE DEVICE



- Do not use the device continuously for more than 30 minutes!
- On first use, it is recommended to start with the minimum intensity for a few minutes. If this test does not generate discomfort, you can adjust it to the maximum intensity.
- In case of discomfort or overheating feelings during use, the device should be turned off. During next use, it is advisable to carry out sessions lasting 10 minutes maximum (3 sessions maximum per day totaling 30 minutes) with a break of 15 minutes between each session.

Power on:

Press  button first to turn it on.

Massage:

Two massage intensities can be selected : low or high intensity.

Note: The last intensity which was operated when last used will turn on by default.

-To get the desired massage intensity, press  button to toggle from low to high intensity.

Heating function:

To activate it, press  button. Press again to turn it off.


**The heat function cannot work alone !
It can only be activated after activating the massage function.**

LED lamp:

To turn it on, press  button. Press again to turn it off.

Note : The LED lamp can be turned on without turning on the device.

POWER OFF:

Press  to turn the device off.

MAINTENANCE AND STORAGE

- Do not immerse the device or any of its parts in water or any other liquids.

Before cleaning the device, make sure that :

- It is turned off,
- It has cooled down if you used before,
- The plug is disconnected.
- Clean the dust out of the device and function remote control with a soft, dry, lint-free cloth.
- Do not use aggressive or abrasive cleaners.
- Cleaning the backrest, cushion and armrests: spot clean with a lightly damp cloth. Then let it dry.
- Only use it when it is completely dry.
- Do not store it in places exposed to direct sunlight, dust or high temperatures.

WARNING: You should not dispose of this device with your household waste. A selective collection system for this type of product is implemented by your local authorities. Please contact your local authorities to find out how and where collection takes place. These restrictions apply because electrical and electronic devices contain dangerous substances that have harmful effects on the environment or on human health and must be recycled.

This symbol indicates that electrical and electronic devices are collected selectively. The symbol shows a waste container crossed out with an X symbol.

Imported by PRODIS 1 rue de Rome 93110 Rosny-sous-Bois/ France/ Made in PRC
Marketed by EUROTUPS VERSAND GMBH D-40764 LANGENFELD

