Gebrauchsanleitung | Instruction leaflet | Instruction d'utilisation  
Instrucciones de empleo | Istruzioni per l'uso | Gebruiksaanwijzing  
Руководство по эксплуатации | Instrukcja użytkowania | Bruksanvisning

# promed NRP-5.4

Nacken-Rückenheizkissen  
Neck-and-back heating pad  
Cousin chauffant pour le dos et la nuque  
Cuscino termico per collo e schiena  
Cojín térmico para las cervicales y la espalda  
Nek-en rugverwarmingskussen  
Электрогрелка для спины и шеи  
Poduszka elektryczna na kark i plecy  
Hals-ryggvärmeryna

DE

EN

FR

IT

ES

NL


RU

PL

SE

AR

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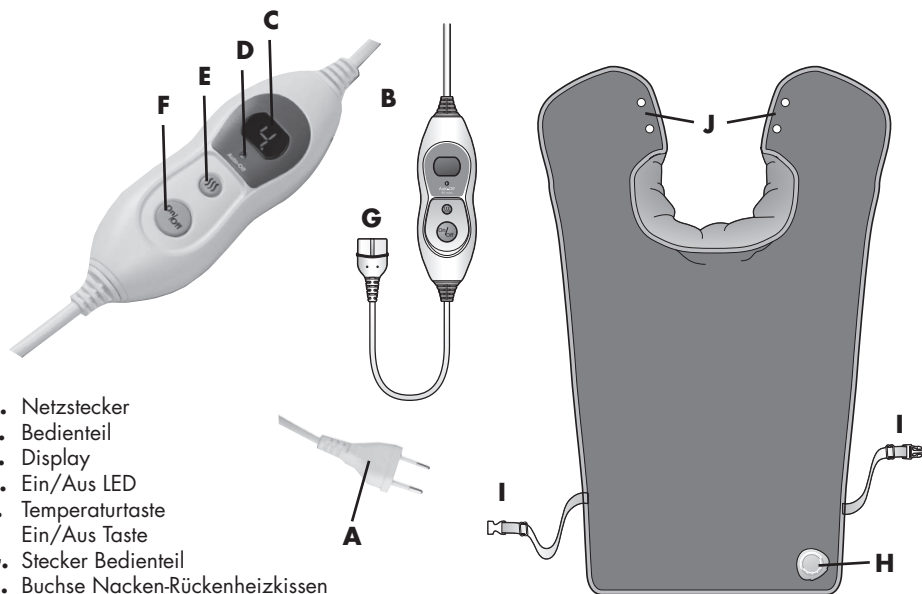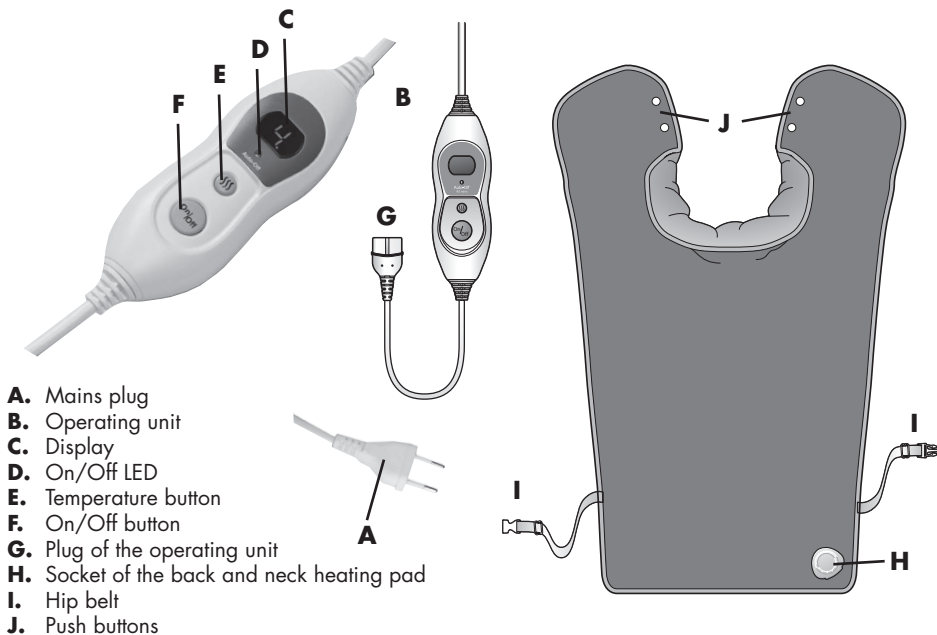
## DESCRIPTION OF THE DEVICE

Please keep this user manual in a safe place!

Please read the instructions for use carefully before using the device for the first time. Thank you for choosing the **Promed NRP-5.4** back and neck heating pad. With the comfortable electronic back and neck heating pad, you can pleasantly warm your body or specific individual regions of your body. The **Promed NRP-5.4** back and neck heating pad adapts to your back and the comfortable belt keeps it in place. Small pockets with ceramic pearls, which are embedded in the heating pad, ensure an even heat distribution. The heating coils integrated in the neck area ensure a completely pleasant sensation of warmth. Thanks to the four temperature settings, you can set the temperature to your individual requirements. The use of this neck and back heating pad promotes the blood flow and thus the relaxation of the tensioned muscle areas and improves the general well-being. You can relax and enjoy its comforting heat whenever you want, for example on your couch. The device is not classified as medical product and thus only intended for domestic use and not suitable for commercial applications or for use in hospitals or doctor's surgeries. To enjoy all benefits of your neck and back heating pad, please thoroughly read these instructions and observe the care and safety instructions. This will enable you to enjoy your heating pad for a longer period of time. Please keep this manual in a safe place for future use or if you want to give the neck and back heating pad to somebody else.

We, the manufacturer, cannot be held liable for any injury or damage to persons or objects resulting from disregard of this manual. We hope you enjoy your new **Promed NRP-5.4**. In the following, we want to improve your understanding of the device. Please read the manual before using the device. This product is equipped with safety functions. Please read the safety instructions thoroughly, however, and only use the product as described in this manual to prevent injury or damage.

DE

EN

FR

IT

ES

NL

RU

PL

SE

AR

# SAFETY INSTRUCTIONS IMPORTANT INSTRUCTIONS.

## RETAIN FOR FUTURE USE.

### Please note the following aspects prior to any commissioning:

- A.** The neck-and-back heating pad may only be used with the controller supplied.
- B.** Voltage supply of the heating pad: 220 - 240 VAC.
- C.** If any problems with the heating pad occur, please return it to the manufacturer.
- D.** Only for indoor use.
- E.** When using an electrical device near children, conscientious supervision is required.

## DANGER!

- Never bring the device into direct contact with water or other liquids.
-   Never use the device in the bathtub, shower, pools, saunas or outdoors.
- Do not touch the device with wet hands.
- Only operate the device with the supply voltage specified.
- The device must not be used outdoors.
- This device uses supply voltage; therefore, the same safety regulations as for any other electrical device apply.
- Prior to any commissioning of the device, check that all components of the device are in perfect condition. In the event of a defect or a malfunction, immediately switch off the device and disconnect the power plug from power supply system.
- Any repair work may only be carried out by the manufacturer. Never open the device. The individual components of this device must not be disassembled, repaired or modified.
- Only use the device in dry rooms and if your body is dry.
- Never operate the device next to other sources of heat.
- Never carry or pull the device by holding it on the power cable. Disconnect the cable from the power socket using the plug and never using the cable.
- Do not pull the power cable over sharp edges or pinch it, do not let it hang down and do not place or suspended it over hot objects (e.g. hobs) or open flames and protect it against heat or oil.
- Only use an extension cable if you previously made sure that it is in a perfect condition.
- If the power cable of the device is damaged, have it immediately replaced by an expert.
- Never use the device after the power cable has been damaged, the device has been dropped or come into contact with water.
- Remove the power plug from the power socket if the device is not in use.
-  The device must not be used by infants and small children (0-3 years).
- This device can be used by children from the age of eight and people with limited physical, sensory or intellectual capabilities or with lack of experience and knowledge if they are supervised by a person responsible for their safety or have received extensive instructions in how to use the device by this person in such a manner that they understand the associated risks.
- Children must not play with the device.
- Cleaning and maintenance work must not be carried out by children without supervision.
-  An electrical device is not a toy. Children do not realise the risks which may result from using electrical devices. Therefore use and store the device out of reach of children.
- Disconnect the device from the power supply system before cleaning it. During cleaning, make sure that no liquids come into contact with the switch and that the device is only used in an absolutely dry condition.
-  Do not insert any needles or pointed objects into the device.
-  READ THE INSTRUCTIONS
- This appliance is not intended for medical use in hospitals.

DE

EN

FR

IT

ES

NL

RU

PL

SE

AR

# SAFETY INSTRUCTIONS IMPORTANT INSTRUCTIONS.

## RETAIN FOR FUTURE USE.

### **WARNING!**

- Never leave the device unattended if children or people unskilled in operating this device are around.
- Only use the device for applications described in these instructions for use.
- Never allow children to play with the packaging material, as there is a risk of suffocation.
- Do not place or use the device on hot surfaces or next to open gas flames to ensure that the housing does not melt.
- Do not store the device next to a washbasin or bathtub, since there is the risk that it can fall or be drawn into the washbasin or bathtub.
- If you suffer from an illness such as varicose veins, thrombosis, diabetes, etc., please consult a doctor before using the device
- Never use the device on open wounds, burns, haematoma, swellings or similar injuries.
- Immediately stop using the device if it feels uncomfortable for you. Heat-sensitive people in particular should use the device carefully and only for short periods of application.
- Under certain circumstances, the electric and magnetic fields emitted by the device may interfere with the function of a cardiac pacemaker. However, these are below the limit values. Please consult your doctor and the manufacturer of your cardiac pacemaker before using the heating pad.  
Electric field strength: max. 5.000 V/m,  
Magnetic field strength: max. 80 A/m  
Magnetic flux density: max. 0.1 mT
- Never use the device on animals.

### **Attention!**

- The device may only be used with the operating unit (HBO5E-4) included in the scope of delivery.
-  Never use the device when it is folded together, creased, or pressed together.
- The device must not be used by small children from the age of three unless the settings have been preset by a parent or legal guardian or unless the child has been sufficiently trained and instructed in how to operate the device.
- Do not use the device when it is wet.
- The device must not be used by helpless and defenceless people or people insensitive to heat.
- Children under the age of three must not use the device due to their inability to react to overheating.
- The electronic components in the switch might heat up when the device is used. Therefore, the switch may not be covered or be placed on the device while it is being operated.
- For storing the device, let it cool down before folding it. Ensure that the device does not crease and that no objects are stored on top of it.

DE

EN

FR

IT

ES

NL

RU

PL

SE

AR

## COMMISSIONING THE DEVICE

Please check that the scope of delivery is complete before commissioning the device.

Ensure that the heating pad is not jammed.

Only use original parts!

Do not use any other operating unit!

### Caution

Please ensure that the heating pad is attached loosely so that some air can circulate between the heating pad and your skin.

Prevent heat accumulation at all costs.

### WARNING – Risk of burning

- A prolonged use of the device at its highest settings can lead to burn injury
- To prevent burn injury to the skin, test the temperature of the device with the back of your hand before using it (do not use the palm of your hand since it is less sensitive than other body parts).
- Do not put the heating pad directly onto your skin. Only use it on your clothing.

### Note:

A slight smell can develop when heating up the device for the first time. This is harmless. If applicable, ensure that it is sufficiently ventilated.

DE

EN

FR

IT

ES

NL

RU

PL

SE

AR

1. Remove the protective film from the display of the control unit.
2. Check if both plugs (plug of the operating unit **(G)** and socket of the neck and back heating pad **(H)**) are connected to each other. If this is not the case, please ensure that they are correctly connected.
3. Put the heating pad on your back and the collar over your shoulders, while the plug of the operating unit **(G)** is directed to the outside. Fix the push buttons **(J)** and the hip belt **(I)** so that the heating pad stays in position.
4. Insert the mains plug **(A)** into the socket. The On/Off LED **(D)** is illuminated and “8” appears on the display **(E)**; both go out after approximately 1 second.
5. Press the On/Off button **(F)** in the display to switch on the neck and back heating pad. “1” appears.
6. Press the temperature button **(E)** to set the desired temperature. You can set the desired temperature in 4 steps). Press the temperature button **(E)** again to increase the desired temperature. Once you have achieved the highest level “4”, the settings start again at “1” by pressing it again. After a few minutes, the heating pad should have reached its maximum temperature. Now you can switch back to a lower heating level.

### Note:

The heating pad is equipped with a safety regulator, which regulates the heat and prevents the device from getting too hot.

7. If you do no longer wish to use the heating pad, press the On/Off button **(F)** to switch off the heating pad and pull the plug out of the socket. Always pull at the plug, not the cable.

### Automatic switch-off

After 90 minutes, the heating pad switches off itself. The ON/Off LED **(D)** is illuminated.

- Activate / deactivate the heating pad with the On/Off button **(F)**.

## TROUBLESHOOTING

Problem	Cause	Error correction
Does not work	The device is not correctly connected to the operating unit	Ensure that the plug of the operating unit <b>(G)</b> and the socket of the heating pad <b>(H)</b> are tightly connected
	The automatic deactivation function has been activated	With the On/Off button <b>(F)</b> you can switch the heating pad on and off.
F is flashing in the display	Did you insert the plug <b>(A)</b> into the wall socket first and then the plug of the operating unit <b>(G)</b> into the socket of the heating pad <b>(H)</b> .	Insert the plug of the operating unit <b>(G)</b> into the socket of the heating pad <b>(H)</b> first at all times before inserting the plug <b>(A)</b> into the socket.
	Has the plug of the operating unit <b>(G)</b> been accidentally removed from the socket of the heating pad <b>(H)</b> .	Pull the plug <b>(A)</b> out of the socket and insert the plug of the operating unit <b>(G)</b> into the socket of the heating pad <b>(H)</b> . Insert the mains plug <b>(A)</b> into a socket again.

## CARE INSTRUCTIONS

### **DANGER - Risk to life due to electric shock**

- Before cleaning the heating pad, pull the plug out of the socket and remove the heating pad from the operating unit.
- Neither the operating unit nor the mains cable may get into contact with humidity, water or other fluids since there is a risk of electric shock. Additionally, the heating pad itself may not get moist while being used. Only use it in a dry state and never use it when it is wet!
-  The heating pad may be washed in the washing machine at 30 °C with mild detergent.
-  Do not iron the heating pad.
-  Do not bleach the heating pad.
-  Do not dry-clean the heating pad.
-  Do not tumble dry the heating pad in a dryer.
- Do not use a high pressure cleaner
- Let the heating pad cool down completely before cleaning it after use.
  1. Pull the plug out of the socket.
  2. Pull the plug of the operating unit out of the socket of the heating pad.
  3. Wash the heating pad at 30 °C in the gentle cleaning cycle. When washing it manually, do not bend or wring out the heating pad.
  4. Let the heating pad air-dry in an unfolded state in a warm environment. Do not use any aids (e.g. clothespins) to put it up for drying. This could damage the internal cables.

### **DANGER - Risk to life due to electric shock**

- After washing, wait until the heating pad has dried completely before inserting the mains plug and connecting it to the power supply. Never use the heating pad in a wet state!
- Never switch on the heating pad for drying.

DE

EN

FR

IT

ES

NL

RU

PL

SE

AR

## CARE INSTRUCTIONS

### CAUTION – Material damage

- Do not dry on a radiator, in the sun or by means of another source of heat.
- Do not iron.
- Do not wring out or deform the heating pad.
- Once the heating pad and the connection socket have dried completely, reinsert the plug of the operating unit into the heating pad socket

### Storage and regular inspection

- Do not crumple the neck-and-back heating pad. Carefully fold the product to store it.
- Store the device in a dry place when it is not in use.
- Allow the device to cool before you fold or roll it up.
- Examine the appliance frequently for signs of wear or damage. If there are such signs, if the appliance has been misused or does not work, return it to the supplier before switching it on again.
- Please do not use any chemical moth protection.
- Do not place any heavy objects on the device during storage, as this may damage the electrical wiring.
- In the event of malfunction or defects, send the device to the manufacturer.

DE

## DISPOSAL

Power tools, accessories, and packaging should be recycled in an environmentally suitable manner.

EN

### Only for EU countries:

Do not dispose of power tools with the domestic waste! In accordance with the European Directive 2012/19/EU on waste electrical and electronic equipment and its implementation into national law, nonfunctioning power tools must be collected separately and recycled in an environmentally suitable manner.

FR

IT

This symbol indicates that disposal of this product in the domestic waste is not permitted within the EU. Waste devices contain valuable recyclable material that should be recycled. Recycling also prevents uncontrolled waste disposal from damaging the environment and human health. Therefore, please dispose of all waste devices in the appropriate collection systems or send the device to the place where you bought the device for disposal. They will then recycle the device.

ES

## WARRANTY

NL

This product was produced with the utmost care and was thoroughly tested before it left the factory. Therefore, we provide a warranty period of 24 months from the date of purchase for this product. **Promed** products conform to their description and the relevant specifications; it is your responsibility to ensure that the product which you have purchased is designed for your intended use.

RU

PL

SE

AR

- In the event of material defects and/or manufacturing errors which appear during proper use and are detected during the warranty period, we will replace all defective parts of the product free of charge within the warranty period, including the labour costs for warranty repairs.
- **The warranty does not cover:**
  - Normal wear on the product
  - Defects which were caused during transportation or storage of the product
  - Defects or damage which were caused through improper use or inadequate maintenance
  - Damage resulting from disregarding the notes in the operating manual
  - Damage resulting from modifications to the product which were not carried out by **Promed**.
  - Damage caused by sharp objects or as a result of torsion, compression, being dropped, an abnormal impact or other actions which are beyond **Promed**'s reasonable control.

## WARRANTY

- Parts which are subject to wear (e.g. moving parts such as ball bearings, etc., fasteners, etc.) are generally excluded from the warranty.
- **A claim under the warranty is not enforceable if:**
  - the product is not sent back in its original packaging or in adequately secure packaging;
  - it has been modified or repaired by a person or company other than **Promed** or one of **Promed**'s authorised dealers;
  - the product has been repaired using parts which are not approved by **Promed**;
  - the serial number / LOT number of the product has been removed, deleted, changed or made illegible.
- For reasons of hygiene, any products which have been exposed to direct contact with a body or with bodily fluids (e.g. blood) must be packaged in an additional plastic bag before they are sent back. In this case, a specific note of this situation must be made in the package or in the accompanying documents.
- For components or products which are replaced during the warranty period, a warranty is only provided for the remainder of the original warranty period; provided that the product or part replacement was carried out by **Promed** or one of **Promed**'s authorised dealers.
- The warranty period begins on the day of purchase. Claims under the warranty must be established within the warranty period. Complaints arising after the expiry of the warranty period cannot be considered.
  - Under these warranty conditions, the warranty shall only apply if the date of purchase can be proven by means of a sales receipt or similar.
  - We reserve the right to make technical and visual changes as well as changes to the features.
- This warranty is only valid and enforceable in the country in which it was purchased by the original buyer, provided that it was **Promed**'s intention that the product should be offered for sale in this country. This warranty is also enforceable in any country within the European Economic Area in which **Promed** has an authorised importer or sales partner. Depending on each country, specific and differing warranties and guarantees may be applicable owing to the particular applicable legislation. These legal provisions are neither excluded nor limited by these warranty conditions. To the extent permissible under national law, the warranty period shall not be extended, renewed or otherwise affected by any subsequent resale, repair or replacement of the product.
- The provisions of the UN Convention on Contracts for the International Sale of Goods shall not apply.
- Our warranty conditions do not affect the seller's legal warranty obligation.
- To the fullest extent possible under the applicable mandatory legislation, these warranty conditions constitute your sole and exclusive legal remedy and shall apply in place of any other express or implied warranty conditions. **Promed** is not liable for any unusual, incidental, punitive or consequential damages including, but not limited to loss of profits, loss of use, loss of revenue, costs of replacement equipment or facilities, third party insurance claims, and damage to property which are the result of the purchase or use of the product or occur owing to a breach of warranty, a breach of contract, negligence, product defects or any other legal or statutory circumstances, even though **Promed** was aware of the possibility of such damages. **Promed** is not liable for any delay in warranty claims.
- **Promed** cannot be held liable for any possible translation errors.

The following information is essential for smooth processing:

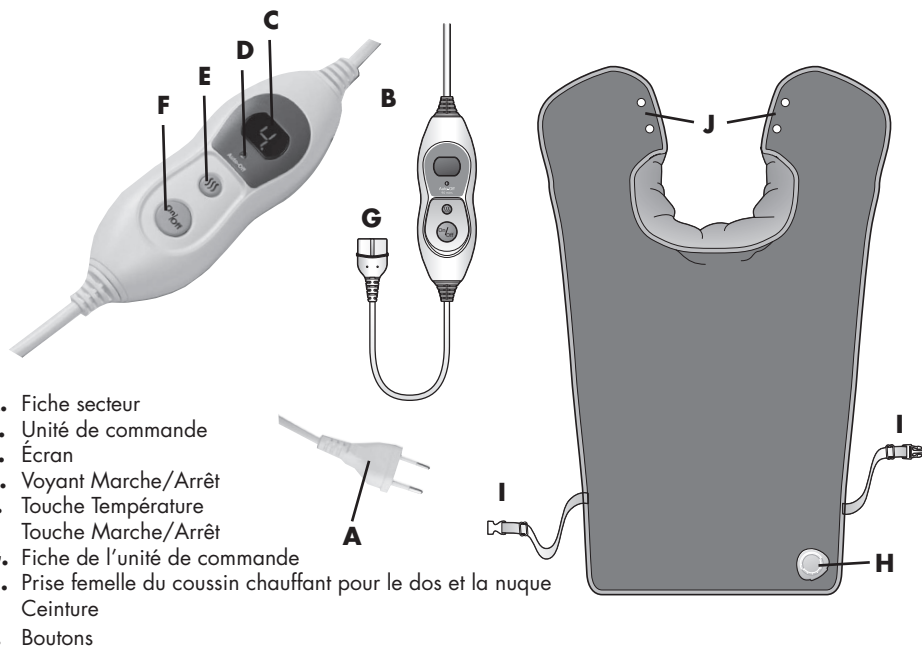
1. Original proof of purchase/receipt or dealer's stamp with date of purchase
2. Defect discovered
3. Product designation / type / serial/batch number

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