
Integrated with both modern technology and traditional medicine, our newly developed and produced cupping instrument is a new health care product, which featured by easy operation, convenient-to-use, low price and is warmly welcomed by suffers.

The traditional cupping treatment has a long history and is widely spread. Cupping is an effective treatment technique rooted from centuries of people battle against diseases. Though there's a folksay that goes "Acupuncture and cupping will cure half of yours disease", the traditional cupping treatment has its disadvantages such as burn, insecurity, limited application scope etc.

Our newly produced cupping instrument replaced the traditional ignition-type air evacuation technique with vacuum air bleed technique, which makes the cupping apparatus absorb on body surface acupoints to cause a contrafluxion or slight blood sludging so as to facilitate smoothness of flows in channels and collaterals, promote blood circulation and treat symptoms such as heatstroke, chill or cold. The cupping instrument can produce a powerful and effective effect in relieving pain, as well as other efficacies such as removing accumulated fat and keeping fit. The cupping instrument can be used for various kinds of treatments.
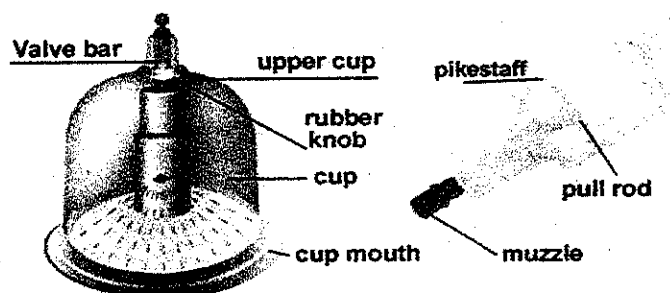
#### 1 Method of Application:

1. pull the upper valve bar of the cup to guarantees smooth air flow before the treatment;
2. fix the vacuum air bleed muzzle to upper part of the cup and then pull the pull rod vertically about 4 times to create appropriate negative pressure;
3. when the treatment has finished, smoothly pull up the valve bar on the upper part of the cup removed the cupping treatment;
4. it is advised to use edible flour dough to make a air tight loop under the cup to ensure air tightness when apply to the part of you with less muscle or other position such as hairy parts, hands and feet as well as arthrosis etc.
5. use disinfectant fluid to wipe the cupping instrument after application. Never soak it.

#### 2 Maintenance and Repairing:

1. avoid crash and scratch to the cup;
2. the rubber knob in the cup shall be kept tightly with the cup. Tighten it with your hand or sharp nose pliers in case a loose rubber knob is detected.
3. you can disassembly the gun block by yourself, and clean out dirt in the gun barrel and then apply some grease in it;
4. you can replace the airtight loop in case of insufficient air bleed before continue to use the product.

#### 3 Methods of Treatment:

1. choose appropriate postures such as face downwards, face upwards, and lie on side of your body according to your state of illness;
2. 10-15 minutes will constitute a treatment, take a treatment everyday or every two days and 7 treatments will constitute a period of treatment. You can start the next period of treatment after a period of treatment. You can start the next period of treatment after a 4-day interval.
3. you can choose from the following treatment methods according to our state of illness and the part of your body: (1) rest cupping (absorb the cup on the body part and let it rest on there of a while); (2) flash cupping (absorb the cup on the body part and then hold the cup to quickly pull it out while a crack voice is heard. This method is only considered to be appropriately adopted when no print mark is caused on the body on the body part.); (3) moving cupping (pierce scarfskin of the focus with a sanitized triangular prism needle to lake it bleed, then pull off the cup at once) etc.
4. the cupping instrument can be used in various forms such as use a single cup of use a number of cups simultaneously, pull off the cup heavily or softly. For elderly, weakling or children, choose to use less cups and pull of the cup more softly and tome of the treatment should be comparatively shorter, generally 5-10 minutes.
5. this product is especially effective for curing heatstroke (i.e. scrape therapy for heat-stroke) and non-fixed acupoints (i.e. any pressure pain point of the body)

#### 4 Precautions:

1. Cupping treatment should not be adopted to serious heart diseases, heart failure, diseases with a tendency of bleeding, patient with tumor, pregnant woman, lower belly of women in menstrual period, ulcer, in case of coma, broken bone or convulsion. Never apply the cups directly on the heart location or on the chest.
  2. The room in which the cupping treatment is done should be kept warm and avoid cold wind. It is advised that you drink some hot water after the cupping treatment and do not wash cupping marks on your body within 30 minutes after the cupping treatment.
  3. If the patient feels swirl, cold sweat, chest distress and other discomfort, removed the cup instantly and drink some warm water. The discomfort will be relieved after a short rest.
- Be sure to apply grease to your skin when adopt moving cupping method. Generally, you can use massage cream.

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