91 x 49

Imported by - Importé par:
GD Import (Genius Ideas & Welllys)
 190 chemin de la Frayère - 06530 Peymeinade - France
 www.gdimport.com
 www.genius-ideas.com
 www.wellys-care.com

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
GB- Training Pad "Next Gen Full Abdo" – PET – 120196

Composition:

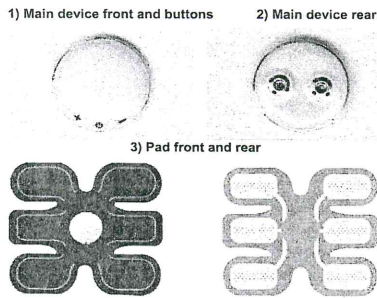Important preliminary information

Refrain from using in the following situations! People who suffer from any of the following diseases or people showing any symptoms should not use the electrostimulator.

- Heart diseases, pacemaker users, persons with implant.
- Abnormal blood pressure, circulatory disorders, varicose veins

- Diabetes, fever, epilepsy, paralysis
- Painful, tuberculosis, nervous diseases
- Menstruation, pregnancy or after giving birth
- Skin diseases, wounds, irritated and/or insensible skin
- Tumors, cancers
- When the body is wet from bathing / showering or sweating
- People receiving medical treatment should consult their doctor prior to using the massager.
- This device is dedicated to be used on belly only. (Abs)
- Do not apply on a painful area.
- Do not apply near to mouth, the head, the neck.
- Do not let the device be in contact with metallic objects such as buckles or jewelry.

CAUTION!

- It is advisable to always start with a low intensity setting with gradual increments to prevent any discomfort. Beginners must stay at a comfortable mode. Caution: Do not assume that by increasing the intensity of stimulation there will be a stronger effect. This may very well not be the case. However, it is possible that excessive stimulation to the skin will result in irritation or red eruptions. Don't use over power, you must always feel comfortable when using this device.
- Be sure not to move the electrode pads without first turning off the power.
- Always keep the electrode pads clean; do not touch the adhesive surface with your fingertips.
- Should you feel any abnormality in either the main unit or your body, stop using immediately.
- Do not use on children, asleep and/or unconscious persons.
- Do not place the electrode pads in the mouth, on or around the genitals, close to the heart, on the chest.
- Do not use the unit in a room with high humidity (bathroom, etc.) or while showering/bathing. Keep a distance least 1 metre to any appliance using short or microwaves.
- Refrain from using when driving. Do not use while sleeping.

To prevent damage

- Do not disassemble the unit, drop it from high places, or subject it to vibrations.
- After use, be sure to turn the power off.
- Cleaning: Separate the pad from the main device after each use, clean the contact surface of the pad with soapy water and let it dry.

Correct use


Installing the battery

- Retire the main device from the pad.
- Gently open the battery compartment behind the front cover of the unit. Helped with your nail or a little flat screwdriver pull on the front cover of the main device in order to open it. Be careful not to force or turn the compartment in order not to break it.
- Insert 1 battery CR2032 type provided and make sure that the (+) is on top of the unit. Then replace the cover until you hear the click sound.

Pad assembly


- Before using it you need to know that the pad must be replaced after 20-30 uses. The loss of adhesiveness, the gel's dryness or dirty grime cannot be removed.
- Connect the main device to the pad as shown on picture.
- After use, attach to the storage sheet.

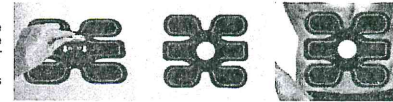
The electro-stimulation pad stimulates the muscles from the abdominal area wherever it is applied. It tones up & increases the muscular volume by weak electrical impulses. This device offers 15 intensities.

How to use it

- Clean the skin where the pad is supposed to be placed. (On abs).
- Place it between your ribs and hips (abdominal area) trying to touch as less as possible the gel.
- And then press it firmly onto your skin.
- Connect the main device to the pad as shown on the picture above.
- Press the button « ON/OFF » to switch on the device. The led light will flash 3 times.

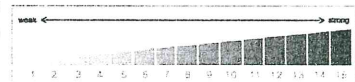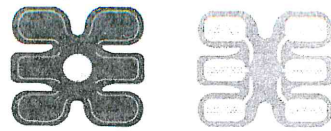6. The device will switch off automatically once the program is completed after 23 minutes. Or if there is an inactivity of 5 sec.
7. To adjust the intensity from level 1 to 15 (1 being the weakest level and 15 the strongest). It is necessary to use the button « + » to increase the intensity and the « - » to decrease.
8. AUTO program: the program works as shown below, during 23 minutes.

- Phase 1: It is the warm up step that lasts 1 minute. The tempo of stimulation increases gradually.
- Phase 2: During 10 minutes the device works and you can choose your comfort intensity.
- Phase 3: It is the rest step that lasts 30 sec. and it allows muscles to rest. (Device is still working).
- Phase 4: It is the fastest step of stimulation that lasts 5 minutes.
- Phase 5: It is the rest step that lasts 30 sec. and it allows muscles to rest. (Device is still working).

- Phase 6: It is the fastest step of stimulation that lasts 5 minutes.
- Phase 7: In this step the tempo slows gradually during 1 minute, and the device turns off at the end of the program.

CAUTION: It is recommended to use this device once a day.

1. Press the button « ON/OFF » to switch of the device. The led light will flash 3 times.
2. Led Indicator:

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Technical information

- Battery: 1 CR2032 battery type provided.
- Frequency 25Hz +/- 10%
- If the device does not work, check if the battery case is firmly closed and the battery is correctly inserted.
- If the intensity is too weak, increase it by pressing the button « + », change the location of the pad on a more muscular abdominal area or simply replace the battery by a new one.
- If the device is faulty, remove the battery for one minute or replace them by new one.

Storage

1. Do not use or store the unit next to magnetic fields or electric waves (TV sets or loudspeakers).
2. Do not place the main unit in areas of high temperature, high humidity, or directly in the sun.
3. Store the unit where there is no moisture.
4. Keep out of reach from children.

Danger Associated with Batteries

- Batteries can be life-threatening when swallowed. Keep batteries and the device out of reach of children. If a battery is swallowed, you must immediately consult medical help.
- Before inserting the battery, check whether the contacts in the device and on the battery are clean and if necessary, clean them.
- Only use the battery type specified in the technical data.
- When changing the battery, check the polarity (+/-).
- Take the batteries out of the unit if they are empty or if you do not use the unit for a prolonged period of time. Thus you will prevent damages that could occur if the batteries leak.
- If a battery has leaked, ensure that the battery acid does not come into contact with your skin, eyes, and mucous membranes. In the event of contact with battery acid, immediately rinse the affected areas off with plenty of clean water and immediately consult a physician.
- Immediately remove a leaked battery from the device. Clean the contacts before you insert a new battery.
- The batteries may not be charged or reactivated with other means, not disassembled, thrown in a fire, submerged in a liquid or short-circuited.

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