Silk'n™

Beautiful Technology

SONICSMILE SG-986

EN USER MANUAL

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

IT MANUALE D'USO

Table of contents

Overview	3
Package contents/device parts	4
General information	5
Reading and storing the operating instructions.....	5
Explanation of symbols	6
Safety	8
Proper use.....	8
Safety instructions.....	9
First use.....	19
Checking the toothbrush and package contents.....	19
Basic cleaning.....	20
Operation.....	21
Charging the rechargeable battery.....	21
Using the toothbrush	22
Brushing your teeth.....	23
Selecting the mode.....	25
Cleaning	27
Storage	28
Troubleshooting.....	29
Technical data.....	31
Declaration of conformity	32
Disposal.....	33
Disposing of the packaging.....	33
Disposing of the toothbrush.....	33

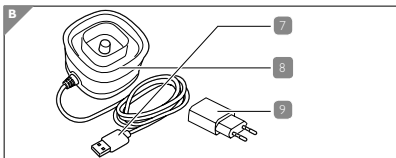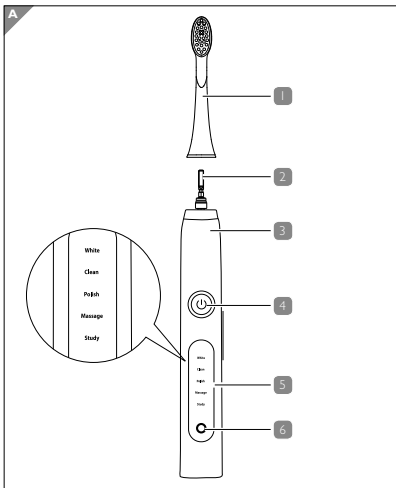Package contents/ device parts

- 1 Brush head, 2×
- 2 Brush attachment point
- 3 Hand piece
- 4 Multifunction button
- 5 LED mode indicator
- 6 Charge level indicator
- 7 Main cord with USB port
- 8 Charging station
- 9 USB adapter

General information

Reading and storing the operating instructions

These operating instructions accompany this SonicSmile SG-986 (referred to below only as the “toothbrush”). It contains important information about start-up and handling.

Before using the toothbrush, please read through the operating instructions carefully. This particularly applies to the safety instructions. Failure to follow these operating instructions may result in severe injury or damage to the toothbrush.

The operating instructions are based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Store the operating instructions for future use. If you pass the toothbrush

on to third parties, please be absolutely sure to include these operating instructions. You can also request these operating instructions in PDF format on our website located at www.silkn.eu.

Explanation of symbols

The following symbols and signal words are used in these operating instructions, on the toothbrush or on the packaging.

This signal symbol/word designates a hazard with moderate degree of risk which may lead to death or severe injury if not avoided.

This signal symbol/word designates a hazard with low risk that, if not avoided, may result in minor or moderate injury.

NOTICE!

This signal word warns against potential damages to property.

This symbol provides you with useful supplementary information on assembly or operation.

Declaration of conformity (see chapter “Declaration of conformity”): Products marked with this symbol meet all applicable Community regulations for the European Economic Area.

Safety

Proper use

The toothbrush is designed exclusively for cleaning teeth. It is only intended for private use and is not suitable for commercial purposes.

Only use the toothbrush as described in these operating instructions. Any other use is deemed improper and may result in damage to property or even personal injury. The toothbrush is not a children's toy.

InnoEssentials accepts no liability for damage caused by improper or incorrect use.

Safety instructions

Risk of electric shock!

A faulty electrical installation or excessive mains voltage may result in an electric shock.

- Only connect the charging station if the mains voltage of the socket corresponds to the specification on the nameplate.
- The main cord of the charging station cannot be replaced. In the event of damage, the charging station must be disposed of.
- Only connect the charging station to an easily accessible socket so that you can quickly disconnect it from the mains in case of a failure.

- Do not operate the toothbrush and the charging station if they exhibit visible damages or the main cord or USB port is defective.
- Do not open the housing; instead, have a qualified professional perform the repairs. Contact a qualified workshop for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection of the device or incorrect operation.
- Only parts that comply with the original device information may be used for repairs. Electrical and mechanical parts, which are essential for providing protection against sources of danger, are

located in this toothbrush.

- The charging station must not be operated with an external timer or with a separate tele-control system.
- Do not immerse the charging station or the main cord in water or other liquids.
- Never touch the USB port with damp hands.
- Never pull the USB port out of the socket by the main cord.
- Never use the main cord as a carrying strap.
- Keep the charging station, toothbrush and USB port away from open flames and hot surfaces.
- Do not kink the main cord and do not lay it over sharp edges.

- Only use the charging station in rooms located indoors.
- Never set up the charging station such that there is a risk of it falling in a tub or sink.
- Never attempt to grab the charging station if it has fallen into water. In such a case, immediately disconnect the USB port.
- Only charge the toothbrush with the enclosed charging station.
- Make sure that children do not insert any objects in the charging station or the hand piece.
- Disconnect the USB port when you are not using the charging station, after the charging process, before you clean the toothbrush or put it away or in the event of

a fault.

WARNING!

Risk of chemical burns!

Fluid leaked from the rechargeable battery may lead to chemical burns.

- Do not allow acid from the rechargeable battery to come into contact with skin, eyes or mucous membranes.
- In the event of contact with acid from the rechargeable battery, immediately flush the applicable areas with plenty of clean water and consult a physician.

WARNING!

Danger of explosion!

Improperly handling rechargeable batteries may cause an explosion.

- Never short circuit the rechargeable batteries.
- Never dismantle rechargeable batteries.
- Never throw rechargeable batteries in fire.
- Never expose rechargeable batteries to excessive heat or direct sunlight.
- Do not in any case open the rechargeable batteries.

**WARNING!**

Danger for children and persons with impaired physical, sensory or mental capacities (e.g. partially disabled persons, older persons with reduced physical and mental capacities) or lack of experience and knowledge (e.g. older children).

- This toothbrush may be used by children and persons with impaired physical, sensory or mental capacities or those lacking experience and knowledge if they are supervised or have been instructed in how to safely use the toothbrush and have understood the risks associated

- with operating it. Children must not play with the toothbrush.
- Cleaning and user maintenance must not be performed by children.
 - Keep children under eight years of age away from the toothbrush.
 - Do not leave the toothbrush unattended while in use.
 - Do not allow children to play with the small parts. Children could swallow small parts and choke.

**CAUTION!****Risk of injury!**

Handling the toothbrush improperly may result in injury.

- Lay the main cord so that it does not pose a tripping hazard.
- Immediately replace worn or damaged brush heads to prevent injury to the mouth and tooth area.

NOTICE!

Risk of damage!

Improper handling of the toothbrush may damage it.

- Do not drop the hand piece. If this does happen, do not use the toothbrush any more. Have it checked by an authorised service point.
- Turn the toothbrush off if you would like to change the brush head.

- Place the charging station on an easily accessible, level and dry surface.
- Ensure that the main cord, the charging station and the toothbrush do not come into contact with hot parts.
- Stop using the toothbrush if its plastic parts exhibit cracks or breaks or are deformed.

First use

Checking the toothbrush and package contents

NOTICE!

Risk of damage!

If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the toothbrush.

- Be very careful when opening the packaging.
 1. Take the toothbrush and its accessories out of the packaging.
 2. Check to make sure that the delivery is complete (see **Fig. A** and **B**).

3. Check whether the toothbrush or the individual parts are damaged. If this is the case, do not use the toothbrush. Contact InnoEssentials at the service address specified on the warranty card.

Basic cleaning

1. Remove the packaging material and all protective foils.
2. Clean all parts of the toothbrush before first use as described in the chapter "Cleaning".

Operation

Charging the rechargeable battery

1. Insert the main cord with USB port **7** into the included USB adapter **9**.
2. Place the toothbrush on the charging station **8**.
The toothbrush will now be inductively charged.
3. The charge level indicator **6** is illuminated in red during the charging process.
4. The toothbrush is fully charged after approx. 16 hours and the charge level indicator goes out.
Once fully charged, the rechargeable battery will keep for up to 3 weeks.

Using the toothbrush

CAUTION!

Risk of injury!

Handling the toothbrush improperly may result in injury.

- Do not apply too much pressure to the toothbrush and your teeth when brushing them.
- Your gums may bleed the first few times you use the toothbrush. If your gums bleed for more than 14 days, please consult a dentist.

- Brush your teeth as close to your gums as possible.
- Brush your teeth using small circular movements.
- Clean up to the spaces between your teeth.

Brushing your teeth

1. Attach a brush head **1** to the hand piece **3** so that the brush attachment point **2** of the hand piece extends into the brush head.
The brush head must audibly and noticeably lock into place.
The brush head must point towards the multifunction button **4**.
2. Hold the brush head in clean water for a moment to wet it.
3. Apply a pea-sized amount of toothpaste to the brush head.
4. Carefully move the brush head up to your teeth.
5. Select the desired mode by pushing the multifunction button within 3 seconds (see chapter "Selecting the mode").
Wait for more than 3 seconds to switch the toothbrush off.

6. Brush the upper and lower parts of your teeth as well as the exterior and interior of your teeth for 30 seconds each without applying pressure.

After every 30 second interval, the toothbrush will stop for 0.5 seconds and automatically switch off after 2 minutes.

7. After each use, clean the toothbrush as described in the chapter “Cleaning”.

Selecting the mode

The LED mode indicator **5** displays the current mode.

Function	Operation	Description
White	Push 1×	For white teeth and less staining.
Clean	Push 2×	For removing plaque, perfect for sensitive teeth and gums.
Polish	Push 3×	For polishing teeth.
Massage	Push 4×	For gently massaging your gums and improving the health of your gums.

Function	Operation	Description
Study	Push 5×	Setting for getting used to the toothbrush. The vibration intensity gradually increases from weak to strong during the first 12 applications.

Cleaning

NOTICE!

Risk of damage!

Handling the toothbrush improperly may result in damage.

- Do not in any case attempt to clean the toothbrush and the charging station in the dishwasher. These parts would be destroyed as a result.
- Do not submerge the charging station in water or other liquids.
- 1. After every use, hold the brush head  under running water and switch the toothbrush on.

2. Remove all leftover tooth paste from the brush head and the hand piece **3** and switch the toothbrush back off.
3. Pull the brush head off of the hand piece.
4. Also clean the brush head by holding it under running water.
5. Dry the brush head and the hand piece.
6. Change the brush head no later than every 3 months.

As soon as the colour of the brush head fades, change it. You can purchase new brush heads on our website at www.silkn.eu.

Storage

All parts must be completely dry before being stored.

- Always store the toothbrush in a dry area.

- Use the original packaging or a similar box for storage.
- Store the toothbrush so that it is securely locked away preventing access by children.

Troubleshooting

Problem	Possible cause	Solution
The toothbrush isn't charging.	The charging station is defective.	Try other sockets or USB connector jacks and charge the toothbrush for at least 16 hours.
	The USB connector jack is defective.	Try other sockets or USB connector jacks and charge the toothbrush for at least 16 hours.
	The charging station is defective.	Replace the charging station within the warranty period.

Problem	Possible cause	Solution
The tooth-brush operates with only a small amount of vibration.	The rechargeable battery is depleted.	Charge the tooth-brush for at least 16 hours.
The tooth-brush does not run for 2 minutes.	The rechargeable battery is depleted.	Charge the tooth-brush for at least 16 hours and make sure that the tooth-brush is properly placed on the charging station.

Technical data

Model:	SG-986
Article number:	SSIPEUW001

Charging station

Mains voltage:	100–240 V~, 50/60 Hz
Power consumption:	1.5 W
Protection class:	III
IP code:	IPX7
Dimensions (L × W × H):	6.2 × 6.2 × 4.3 cm

Hand piece

Rechargeable battery:	3.7 V, 700 mAh
Operating time:	42 times/full charge
Timer:	0.5 second pause after approx. 30 seconds, deactivation after 2 minutes

Dimensions

(L × W × H): 2.6 × 2.7 × 25.4 cm

Adapter

Model number: ZD5C050060EUDU

Input: 100–240 V~,
50/60 Hz, 0.2A

Output: 5.0V --- 600mA

Protection class: II

Declaration of conformity

The EU declaration of conformity can be requested from the address specified on the warranty card.

Disposal

Disposing of the packaging

Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

Disposing of the toothbrush

(Applicable in the European Union and other European states with systems for the separate collection of reusable waste materials)

Old devices must not be disposed of with household waste!

If at some point it is no longer possible to use the toothbrush, each consumer **is required by law to dispose of old devices separately from**

their household waste, e.g. at a collection point in their community/borough. This ensures that old devices are recycled in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries do not belong in the household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or borough or to a retailer, so that they can be disposed of in an environmentally friendly manner.

Bring the entire toothbrush (with the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

*labelled with: Cd = cadmium, Hg = mercury, Pb = lead

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InnoEssentials Deutschland GmbH