prorelax[®] Therapy

15100

ThermoTENS

DE GEBRAUCHSANLEITUNG
UK OPERATING INSTRUCTION
FR MODE D'EMPLOI

prorelax[®] Therapy

ThermoTENS

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ENGLISH

IMPORTANT INSTRUCTIONS! KEEP FOR LATER USE!

Read the instructions for use carefully, especially the safety instructions, before using the device and keep the instructions for use for future use. If you pass the device on to third parties, be sure to include these instructions for use.

TABLE OF CONTENT

Warning Signs, Symbols and their Meaning	17
Intended Use	18
Important Safety Instructions	18
Parts and Components	22
Delivery Content	22
Before first use	23
Usage	24
Cleaning and Care	26
Technical Data	26
Recycling	27
Guarantee Provisions	27

Warning Signs, Symbols and their Meaning

	These operating instructions belong to this device. It contains important information on commissioning and handling. Read this instruction manual completely. Failure to follow this instruction could result in serious injury or equipment damage.
	WARNING: Failure to heed the warnings could result in injury to the user.
	DANGER: If the instructions are ignored, there is a risk of damage to the device.
	NOTE: Observe these notes for preparation or operation.
	Do not use if you have implanted electrical devices (e.g. pacemakers).
	Do not use near the heart.
 Operating	Permissible temperature and humidity during operation.
 Storage/Transport	Permissible temperature and humidity during storage and transport.
	Expiry Date
	Dispose of the packaging in an environmentally friendly manner.
	This device can provide RMS outputs exceeding 10 mA, averaging over each 5 second interval.
	CE - Certified according to MDD Directive 93/42/EC Medical Devices
	Device classification: Type BF
	Manufacturer
	Manufacturing date
	LOT number
	Disposal in accordance with the Waste Electrical and Electronic Equipment (WEEE) and Battery Directive.

Intended Use

The device is designed exclusively for use in private households and may not be used for commercial or medical purposes.

We assume no liability for damage to the product, damage to property or personal injury resulting from careless, improper, incorrect use of the product or use not in accordance with the purpose specified by the manufacturer. This device is intended for use in adult to temporarily relieve back pain, joint pain, neuralgia, pain after injuries to musculoskeletal system, pain with circulatory problems, and chronic pain through various causes.

Important Safety Instructions

Use of the device does not replace medical advice or treatment.

So always consult your doctor first.

WARNING!

In order to prevent damage to health, no electrical stimulation device shall not be used in the following cases or on the named parts of the body:

- With implanted electrical devices (e.g. cardiac pacemakers)
- In the presence of metal implants
- For insulin pump wearers
- If you have a high fever (e.g. $> 39\text{ }^{\circ}\text{C}$)
- In the case of known or acute cardiac arrhythmias and other cardiac excitation and conduction disorders
- For seizure disorders (e.g. epilepsy) If you are pregnant
- In the presence of cancer
- After surgery where increased muscle contractions could disrupt the healing process
- An application near the heart must be avoided. Stimulation electrodes should not be applied to any part of the anterior chest (indicated by the ribs and breastbone), especially not on both of the large pectoral muscles. Here it can increase the risk of ventricular fibrillation and cause cardiac arrest.
- On the cranium, in the area of the mouth, pharynx or larynx

- In the neck / carotid artery area
- In the genital area
- On acutely or chronically diseased (injured or inflamed) skin (e.g. in the case of painful and painless inflammation, redness, skin rashes (e.g. allergies), burns, bruises, swelling and open wounds as well as those in the healing process, on surgical scars which are in the process of healing are understood)
- In high humidity environments such as in the bathroom or while bathing or showering
- Do not use after drinking alcohol
- With a simultaneous connection to a high-frequency surgical device.
- Acute or chronic diseases of the gastrointestinal tract.
- Stimulation should not be applied over or through the head, directly on the eyes, at the mouth, on the front of the neck (especially not to the carotid artery), or with pads placed on the chest and upper back or across the heart.

Before using the device, consult your doctor if any of the following applies to you:

- Serious illnesses, in particular if you suspect or have been diagnosed with high blood pressure, a blood coagulation disorder, propensity to thrombo-embolic conditions or recurrent malignant growths
- Any skin conditions
- Unexplained chronic pain in any part of the body Diabetes
- Any sensory impairment that reduces the feeling of pain (e.g. metabolic disorders)
- If you are receiving medical treatment
- In the event of complaints linked to stimulation treatment
- If you suffer from persistently irritated skin due to long-term stimulation at the same electro-
de site.

WARNING!

Only use digital TENS units:

- On adults.
- For the intended purpose and as specified in these instructions for use. Improper use can be dangerous.
- For external use.
- With the original accessories supplied, which can be re-ordered. Failure to do so invalidates the warranty.

PRECAUTIONS:

- Always pull firmly on pads to remove them from the skin to prevent injuries in the unusual case of highly sensitive skin.
- Keep the device away from sources of heat and do not use it in close proximity (approx. 1 m) to shortwave or microwave devices (e.g. mobile phones), as doing so can result in unpleasant current peaks.
- Do not expose the device to direct sunlight or high temperatures.
- Protect the device from dust, dirt and humidity.
- Never immerse the device in water or other liquids.
- The device is suitable for self-treatment.
- For hygiene reasons, the electrodes may only be used on one person.
- If the device does not work properly, or you feel unwell or experience pain, stop using it immediately.
- Switch off the unit or the respective channel first before removing or relocating electrodes to prevent unintentional stimulation.
- Do not modify electrodes (e.g. by cutting them). This increases the current intensity, which is potentially hazardous (max. recommended output value for the electrodes is 9 mA/cm², an effective current intensity beyond 2 mA/cm² requires increased awareness).
- Do not use the device whilst asleep, driving a vehicle or operating machinery.
- Do not use whilst undertaking any activity where an unexpected reaction (e.g. strong muscle contractions even at low intensity) could be dangerous.
- Ensure that no metallic objects (e.g. belt buckles or necklaces) come into contact with the electrodes during stimulation. If you are wearing jewellery or have piercings in the area to be treated (e.g. a navel piercing), these must be removed before using the device. Failure to do so could result in spot burns.
- Keep the device away from children.
- Make sure not to confuse the electrode cables including contacts with your headphones or other devices and do not connect the electrodes to other devices.
- Do not use the device whilst using other devices that transmit electrical impulses into your body.
- Do not use in the vicinity of highly flammable substances, gases or explosives.
- Do not use rechargeable batteries and always use the same battery types.
- During the initial few minutes, use the device while sitting or lying down to minimise the risk of injuries as a consequence of isolated cases of vagal responses (feeling of faintness). If you feel faint, immediately switch off the device, lie down and support the legs in an elevated position (approx. 5 – 10 min).
- Treatment of the skin with moisturising lotions or ointments beforehand is not recommended as this considerably increases the electrode wear and may cause unpleasant current peaks.
- This device is not intended for use by children or people with restricted physical, sensory (e.g.

reduced sensitivity to pain) or mental skills or a lack of experience and/or lack of knowledge, unless they are supervised by a person who is responsible for their safety or are instructed by such a person in how to use the device.

- If the adhesive capability of the adhesive pads decreases, please replace them immediately. Only use the device again with new adhesive pads. Otherwise the unequal adhesion of the adhesive pads may lead to skin injuries. Replace the pads with new ones after having used them 20 times at the latest.

Damage

- Check the device for signs of wear and tear or damage. If there are such signs of wear and tear or damage or if the device was used improperly, it must be returned to the manufacturer or retailer before further use.
- Switch the device off immediately if it is faulty or not working properly.
- Do not attempt to open and/or repair the device yourself. Repairs may only be carried out by Customer Services or authorised retailers. Failure to comply will result in voiding of the warranty.
- The manufacturer is not liable for damage resulting from improper or careless use.

Battery safety instructions


- **WARNING:** To charge the battery, only use a suitable USB charger that meets the requirements of the specific standards.
- Never expose batteries to excessive heat (e.g. direct sunlight, fire) and never throw them into fire. The batteries could explode.
- The batteries are built into the product and cannot be replaced. Before disposing of the product, contact an electrician or your dealer to remove the batteries and dispose of them in an environmentally friendly manner.
- Never damage the batteries. There is a risk of explosion and fire if the battery case is damaged! Never short-circuit the contacts/terminals of the battery. Do not throw the battery or the product into a fire. There is a risk of explosion and fire!

Note:

This device can be used to reduce pain (**TENS = transcutaneous electrical nerve stimulation**). TENS means electrical nerve stimulation through the skin. Chronic pain can be treated with different methods, e.g. with medication, surgeries but also with other methods without side effects such as acupuncture or TENS.

Parts and Components

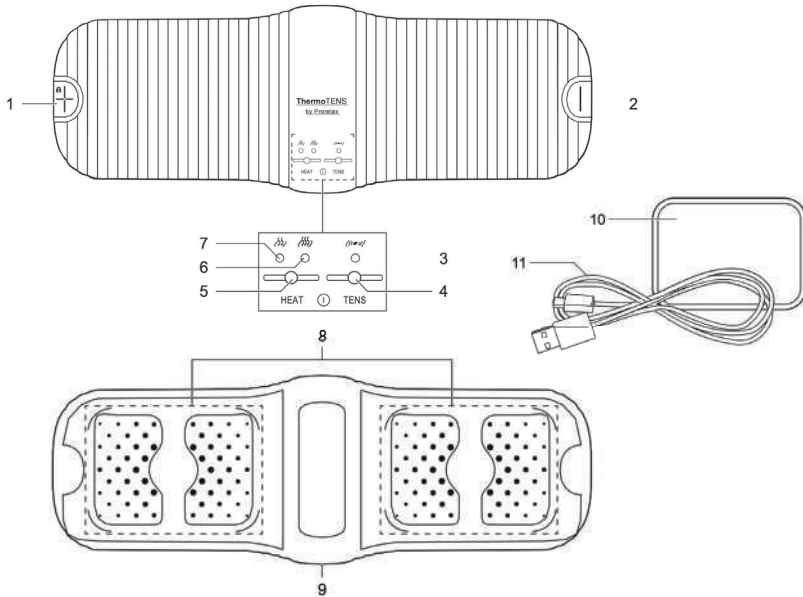- 1. Intensity button **+**, button lock
- 2. Intensity Button **-**
- 3. TENS indicator LED
- 4. Function/pause button **TENS**
- 5. Heating function button **HEAT**
- 6. Heat Level Indicator LED, High
- 7. Heat Level Indicator LED, Low
- 8. Electrodes
- 9. USB C charging socket
- 10. Gel sticky pads (2x)
- 11. USB cable

Delivery Content

1 x TENS device
 1 x USB charging cable
 2 x Gel sticky pads

Before first use

Charging the battery

Allow the device to charge for at least 2 hours before using it for the first time.

1. Connect the supplied USB cable to a USB power source and to the USB C socket of the device.
2. Plug the USB charger (5V/1A) into a suitable socket.
3. During the charging process, the TENS indicator LED lights up green. When the battery is fully charged, the TENS indicator LED turns off.

You cannot use this device while it is charging. Only use this product with a suitable USB cable and a suitable USB power source that meets the specific applicable standards.

Attaching the gel pads

- Remove one of the protective films with care (Fig. 1).
- Attach the gel pads to the electrode area and ensure that the edge of the gel pad does not protrude outside the electrode profile. Only attach one gel pad to each side so does it covers both electrodes on each side (Fig. 2).
- Remove the second protective film slowly and carefully to make sure the gel pads are not damaged.

Use this product only with the original gel pads supplied, which can be re-ordered. Failure to do this invalidates the warranty. The manufacturer is not liable for damage resulting from improper or incorrect use.

Fig.1


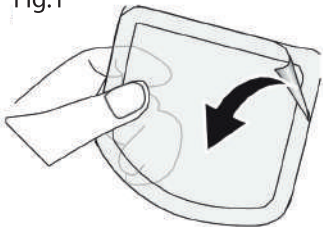
Fig.2

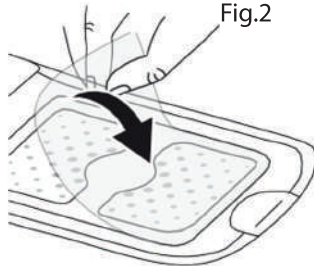Usage

- There are total 3 modes (TENS mode, TENS and heating mode combined, heating mode) available. This device has two levels of heat function, and 15 intensity levels for TENS application.
 - We recommended you start using a low intensity levels until you have got used to the device.
 - This device has an automatic switch-off function. It switches off automatically after 30 minutes of using combined TENS and heat mode or after 20 minutes of using heat only or TENS only.
- ① Please stick the device on the painful area with gel pad. This device has a function which detects whether the device is in contact with the body properly, if there is not sufficient contact with the body, the device beeps and switches itself off.

② **Choose one of the following modes:**

A) TENS-Mode:

- Press the TENS button or the **+** button for **1 second** to turn on the device and start the TENS function. **The TENS indicator lights up green.**
- Briefly press the **+** button to increase the intensity. Briefly press the **-** button on the device to decrease the intensity.
- **Keylock:** Press and hold the **+** button for **2 seconds** to activate or deactivate the key lock. An acoustic signal beeps when the key lock is activated and deactivated.
- If the key lock is activated, the TENS intensity or the heating level cannot be adjusted, but the device can be switched off by pressing and holding the **-** button.
- If necessary, briefly press the TENS button to pause the TENS mode or to resume it.
- If necessary, briefly press the HEAT button to switch to the combined TENS and heating mode (TENS and heating indicator LEDs light up).

B) Heating mode

- **Long press the HEAT button** to start the low heat setting. The corresponding indicator LED lights up green.
- Press the HEAT button again to get to the high heating level. The corresponding indicator LED lights up green.
- **Keylock:** Press and hold the **+** button for **2 seconds** to activate or deactivate the key lock. An acoustic signal beeps when the key lock is activated and deactivated. If the key lock is activated, the TENS intensity or the heating level cannot be adjusted, but the device can be switched off by pressing and holding the **-** button.

- If necessary, briefly press the TENS button to switch to the combined TENS and heating mode (TENS and heating indicator LEDs light up).

C) Combined TENS and heating mode

- Press the button **TENS and HEAT together** to activate the combine mode. Alternative you can first press the button TENS and then the HEAT button. The corresponding indicator LED lights up green.
- Briefly press the **+** button to increase the intensity. Briefly press the **-** button on the device to decrease the intensity.
- Keylock: Press and hold the **+** button for 2 seconds to activate or deactivate the key lock. An acoustic signal beeps when the key lock is activated and deactivated.
- If the key lock is activated, the TENS intensity or the heating level cannot be adjusted, but the device can be switched off by pressing and holding the **-** button.
- If necessary, briefly press the TENS button to pause the TENS mode or to resume it.
- If necessary, briefly press the HEAT button again to switch from low heat level to high heat level. If the heating function is activated, a long press on the HEAT button turns off the heating function.

- ③ As previously described, the device will turn off automatically after either 20 or 30 minutes of operation. The device can be switched off manually by pressing and holding the **-** button or the TENS button.

Attention!

This device has a function which detects whether the device is in contact with the body properly, if there is not sufficient contact with the body, the device beeps and switches itself off. If the device is attached to your body but still beeps, clean the gel pads and try starting the device again. If the continuous beeping persists, replace the gel pads.

Notes on use

- Switch off the device before removing it.
- Wait at least 30 minutes between applications.
- Only use the device as often when you find comfortable. If treatment with this device starts to feel unpleasant, stop using it.
- If the adhesive capability of the gel pads decreases, replace them immediately. Do not use the device again until the gel pads have been replaced, as the uneven adhesion of the gel pads could lead to skin lesions.
- In the event of unpleasant current peaks, switch off the device and move it slightly.

Cleaning and Care

Danger!

Never use solvents or scouring agents, hard brushes, metal or sharp objects for cleaning. Solvents are harmful to health and attack the plastic parts; abrasive cleaners and auxiliaries will scratch the surface.

- Clean the outer surfaces of the product with a slightly damp cloth and then dry them thoroughly. Gently clean the adhesive gel pads with a damp, lint-free cloth.
- Store the product in a dry, cool place away from moisture and direct sunlight.
- Store the product in a place inaccessible to children.

Technical Data

Type and model	Type: FT-810 / Model: 15100
Output waveform	Biphasic rectangular pulses
Pulse length	200 - 400 μ s
Pulse frequency	6 - 45 Hz
Output voltage	Max. 80 V _{p-p} (at 1 k Ω)
Output current	80 mA _{p-p} (at 1 k Ω)
Heat mode settings	HI: up to 43 °C LOW: up to 41 °C
Rated voltage	3,7 V ---
Rate input voltage / current	5 V --- , 1 A
Rechargeable battery	Li-Ion, 3.7 V --- , 1000 mAh
Treatment time, automatic power off	TENS mode: 20 minutes, TENS and heat mode: 30 minutes, Heat mode: 20 minutes
Intensity settings	0 to 15
Operating conditions	5-40 °C at a relative humidity of 15-90%, atmospheric pressure: 700 hPa to 1060 hPa
Storage conditions	0-40 °C at a relative humidity of 0-90%
Dimensions	100 mm \times 300 mm \times 15 mm
Weight	102 g
Gel pad life cycles	20

Recycling

This marking means that this product should not be disposed of with your general household waste. To prevent harm to human health and the environment due to uncontrolled waste disposal, recycle the product responsibly, promoting sustainable reuse of material resources. To return your old device, please use the return and collection systems or contact the point of sale where you purchased the product. There, this product is sent for environmentally friendly recycling. You are legally obliged to return old electronic devices after use. This disposal is free of charge for you. Small old electrical appliances (edge length up to 25 cm) can be purchased from dealers with a sales area for electrical and electronic equipment of at least 400 m² or Food retailers with a total sales area of at least 800 m² that offer electrical and electronic equipment at least several times a year can be returned free of charge (without buying a new device). You are obliged to take back old

electrical and electronic equipment free of charge. When buying a new device of the same type of device, which essentially fulfills the same functions as the new device, larger old devices can be purchased at a corresponding dealer can be returned free of charge. For the conditions of returning an old device in case of delivery of the new device, please contact your dealer or point of sale. In the case of sales using long-distance means of communication, all storage and shipping areas are considered sales areas of the distributor. Before disposing of the product, remove all batteries and rechargeable batteries and all lamps (if included) that can be removed non-destructively.

We would like to point out that you are responsible for deleting personal data (if included) on the product to be disposed of.

The symbol of the crossed-out wheeled bin on batteries or accumulators means that they must not be disposed of with household waste at the end of their service life. If batteries or accumula-

tors contain mercury (Hg), cadmium (Cd) or lead (Pb), you will find the respective chemical symbol below the crossed-out garbage can symbol. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a retail store or at another collection point near you. Addresses of suitable collection points can be obtained from your city or local government. Batteries may contain substances that are harmful to the environment and human health. Particular caution is required due to the special risks involved in handling lithium-containing batteries. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.

Please avoid generating waste from old batteries as much as possible, e.g. by giving preference to batteries with a longer service life or rechargeable batteries. Please avoid littering public spaces by not leaving batteries or battery-containing electrical and electronic devices lying around.

Guarantee Provisions

EUROMEDICS GmbH, Hauptstraße 169, 59846 Sundern, Germany (hereinafter called "EUROMEDICS") provides a warranty for the prorelax® ThermoTENS (hereinafter called the "Product") subject to the following preconditions and in the scope described below.

The warranty conditions described below do not affect the legally prescribed warranty obligations of the seller from the contract of purchase with the buyer. Moreover, the warranty applies irrespective of mandatory legal liability provisions, such as those according to the product liability law, or in cases of mala fide intent or gross negligence, or owing to injury to life and limb, or to health by EUROMEDICS or its agents.

EUROMEDICS guarantees the defect-free operability and the completeness of the product, especially that the product is free of material defects, manufacturing defects and construction defects. The decisive factor in this regard is the state of the art of the science and technology at the time of manufacturing. The product must have manifested the defect that caused the damage at that time itself. The warranty period is 2 years from the date of purchase of a new, unused product by the purchaser.

The warranty applies worldwide, i.e. regardless of the country in which EUROMEDICS has marketed the product and the purchaser has purchased it.

This warranty only applies to products that the purchaser has procured as a consumer and exclusively for personal purposes in the framework of domestic use (i.e. not for products that were procured or used in the framework of a commercial or independent professional activity).

If the product is found to be deficient in operability in accordance with the following provisions during the warranty period, EUROMEDICS will provide a free replacement according to these warranty conditions.

If the purchaser wishes to file a warranty claim, he must first contact EUROMEDICS customer service:

EUROMEDICS KUNDENSERVICE
Biberweg 12
D-53842 Troisdorf
Tel: +49 (0)2208/9217999
Germany
E-mail: info@euromedics.de

The purchaser will then get further information on how to process the guarantee case, for example, where he can send the device free of cost and what documentation is required.

Filing a warranty claim can only come into consideration if the buyer can submit

- a copy of the invoice/proof of purchase and
- the original product

to EUROMEDICS or to a service centre nominated by EUROMEDICS upon demand.

Expressly excepted from this warranty are

- products that are used, cleaned, stored or maintained incorrectly and/or not in compliance with the conditions laid down in the operating manual, as well as products that were opened, repaired or modified by the purchaser or by a service centre not authorised by EUROMEDICS;
- products that were sold as second quality or used articles;
- consequential damages arising from a defect in the product (however, in this case, there may be entitlements from product liability or from other mandatory legal liability regulations).

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Type and model	Type: FT-810 / Model: 15100
Output waveform	Biphasic rectangular pulses
Pulse length	200 - 400 μ s
Pulse frequency	6 - 45 Hz
Output voltage	Max. 80 V _{p-p} (at 1 k Ω)
Output current	80 mA _{p-p} (at 1k Ω)
Heat mode settings	HI: up to 43 °C LOW: up to 41 °C
Rated voltage	3,7 V ---
Rate input voltage / current	5 V --- , 1 A
Rechargeable battery	Li-Ion, 3.7 V --- , 1000 mAh
Treatment time, automatic power off	TENS mode: 20 minutes, TENS and heat mode: 30 minutes, Heat mode: 20 minutes
Intensity settings	0 to 15
Operating conditions	5-40 °C at a relative humidity of 15-90%, atmospheric pressure: 700 hPa to 1060 hPa
Storage conditions	0-40 °C at a relative humidity of 0-90%
Dimensions	100 mm \times 300 mm \times 15 mm
Weight	102 g
Gel pad life cycles	20

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HIVOX BIOTEK B.V.